

VDW.GOLD®

Multifunctional endo motor with integrated apex locator



Versatile in function and range of application

- Choice of file systems
- Pre-programmed torque and speed settings for all major NiTi systems
- Great flexibility due to easy storage of individual torque and speed settings
- Intelligent control unit with automatic safety functions
- Separate length determination with manual instrument or simultaneous length control during preparation
- Convenient battery operation; motor can also be used while charging


The VDW logo consists of three blue diagonal lines to the left of the letters 'VDW' in a bold, blue, sans-serif font.

Endo Easy Efficient®



Functional control panel

- ANA** increases safety for difficult anatomy by reducing torque limit
- ASR** automatic STOP and reverse rotation when torque limit is reached
- MENU** simple and practical menu structure
- CAL** calibration of newly sterilised contra-angle

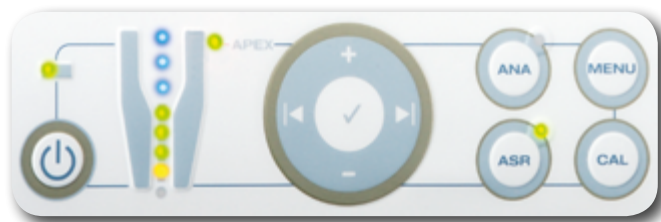

Integrated apex locator

- Simultaneous length control during preparation or separate length determination with manual instrument
- Auto-STOP upon reaching apex can be activated
- Continuous monitoring of the file position via LED display
- Reliable length determination without X-ray



High quality contra-angle with integrated apex locator contact

- Small head for good visibility and easy handling
- No complicated file contact clip
- Fast file change



Range of products

- Ref. 1103 **VDW. GOLD®**
Including contra-angle and foot control
- Ref. 1079 **Additional 6:1 contra-angle with small head**
- Ref. 1105 **VDW. GOLD® (Ref. 1103) + Mtwo® StarterKit**
Contents Mtwo® StarterKit:
2 Mtwo® of sizes 10/.04, 15/.05, 20/.06, 25/.06 each
1 Mtwo® of sizes 30/.05, 35/.04, 40/.04, 25/.07 each
1 Mtwo® SystemBox
3 plastic training blocks
2 x 3ml FileCare®EDTA
1 user CD-ROM and instruction brochure
- Ref. 1104 **VDW. GOLD® (Ref. 1103) + FlexMaster® StarterKit**
Content FlexMaster® StarterKit:
Each 1 FlexMaster® BasicKit 21 and 25mm
1 FlexMaster® SystemBox
3 plastic training blocks
2 x 3ml FileCare®EDTA
1 instruction brochure